

# StatSpin Express 3 Centrifuge



## Faster turnaround times and greater productivity

With the StatSpin Express 3, your laboratory will experience faster turnaround times and greater productivity. This primary tube centrifuge holds up to 8 tubes per run, with the ability to run 10 mL tubes or 5 mL with an adapter.

The StatSpin Express 3 is simple to use and patient samples can be run on demand without batching runs. The StatSpin Express 3 features three settings that provide results in 2, 3, or 5 minutes without the need for speed adjustment.

### Key Product Features:

#### Safety Features

- Lid remains locked until the rotor comes to a complete stop.

#### Compact

- Small, round footprint takes up little bench space and fits well in small corners.

#### Ease of Use

- Just load tubes, select time, and press start.




## StatSpin Express 3 Centrifuge Specifications<sup>1</sup>

<b>Dimensions:</b>	Height 7 in (18 cm) Width 11 in (26.7 cm) Depth 10 in (25.4 cm) Weight 8.75 lbs (4.3 kg)
<b>Cycle Times/Speed/RCF:</b>	2 minutes / 7200 rpm / 4400 x g (+/-5%) 3 minutes / 7200 rpm / 4400 x g (+/-5%) 5 minutes / 5600 rpm / 2685 x g (+/-5%)
<b>Acceleration Time:</b>	≤ 25 seconds to 95% of speed
<b>Deceleration Time:</b>	≥ 30 seconds
<b>Tube Sizes:</b>	Accepts blood tubes 1.5-10 mL
<b>Power Requirements:</b>	100-240 V, 50/60 Hz, 3 amp
<b>Warranty:</b>	Two years

### Ordering Information

SSX3

StatSpin Express 3 Centrifuge for 100-240 V, 50/60 Hz.  
Includes 8-place rotor, removable bowl liner, eight 5 mL inserts  
and power cord.

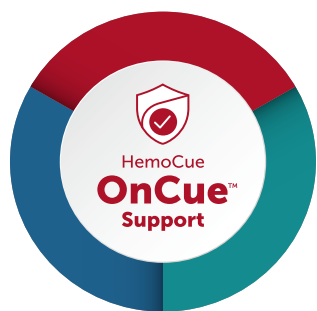

### HemoCue OnCue™ Support

HemoCue OnCue™ Support is a comprehensive program designed to meet your needs as quickly and conveniently as possible.

Powered by HemoCue employees, you get answers and education from the experts who know our products best.

At your office, online, or on the phone, we're ready to do more for you.

<sup>1</sup>StatSpin Express 3 Centrifuge Instructions For Use 55-005187-001.EC

Specifications subject to change without prior notice. The StatSpin Express drive and suspension system are covered under U.S. Patent No. 5,566,919.

Certified by a Nationally Recognized Testing Laboratory (NRTL) for the US and Canada.

©2020 Beckman Coulter. All rights reserved.

StatSpin products are developed and manufactured by Beckman Coulter, Inc. and sold through HemoCue America, a division of Radiometer America Inc. Beckman Coulter, the stylized logo, and the Beckman Coulter product and service marks mentioned herein are trademarks or registered trademarks of Beckman Coulter, Inc. in the United States and other countries.

HemoCue America | California, USA | Phone 800.881.1611 | Fax 800.333.7043

web@hemocue.com [hemocue.us](http://hemocue.us)

©2020 HemoCue America. All rights reserved. MMUS-0054[3]

